


Included in the Box

① Neuvago Wellness Device




② Earpiece (Left) & Ear Clip

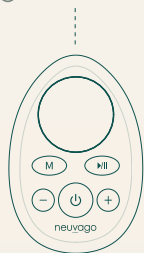

③ USB-C Charging Cable



④ User Manual & Quick Start Guide

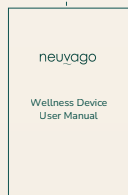

⑤ Conductive Gel Tube



⑥ Spray Bottle



neuvago

Quick Start Guide



1

Connect and power on the device

Plug the earpiece cable into the top port, then press and hold the power button for 3 seconds to turn on.



Power button

2

Prepare and position the electrodes

Moisten the electrodes with clean water or gel.
Start with clean skin and adjust as needed for your comfort.



3


Properly fit the earpiece

Place the stimulation tips on the inner areas shown below on your LEFT ear and secure comfortably.



Correct placement

Wrong positioning may bring uncomfortable feelings or pain.

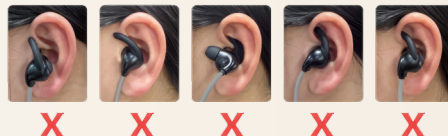

Placement guidance

The unit will show a green symbol of the earpiece when the electrodes are correctly placed. The color will be orange if you have poor contact or disconnected. Ensure the electrode is moistened with water or gel for contact.



4

Adjust the intensity

Press + or - to find your ideal level of comfort.



Scan QR-code for instructions

NOTICE

This device is designed for general wellness and relaxation.
If the sensation feels too strong or uncomfortable, lower the intensity or pause your session and adjust as needed.
Use on clean, dry skin for the best experience.
If you have any questions, please see our website.