

Neuvago Wellness Device

User Manual





June 2026

<https://www.neuvago.com>

Neuvago Device Pack

This pack contains:


- 1 x Neuvago Wellness Device
- 1 x Earpiece (Left)
- 1 x USB-C Charging Cable
- 1 x Felt Case
- 1 x User Manual
- 1 x Quick Start Guide
- 1 x Spray Bottle
- 1 x Gel Tube

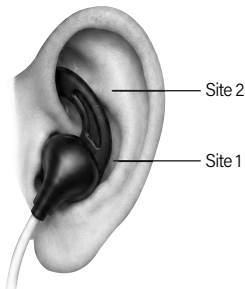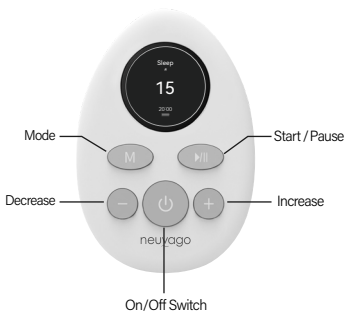

Precautions

- The earpiece must be in contact with the ear.
- Avoid using the earpiece on wounds, scars, bruises, abrasions, or inflamed skin.
- Do not place the earpiece in front of or behind the heart during use.

Device Overview

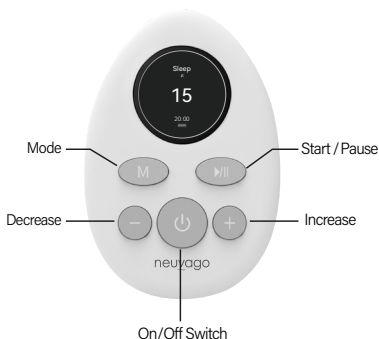

Button functions

- **Mode (M):** Switch between device modes (Sleep, Relax, Meditation and Relief).
- **Start / Pause (▶||):** Start, pause, or resume a session.
- **Increase (+):** Increase the intensity level .
- **Decrease (-):** Decrease the intensity level.
- **Power(⏻):** Turn the device on or off.

Before your first session

1. Charge the device fully before use.
2. Connect the earpiece to the device.
3. Moisten the earpiece contact area with a few drops of water.
4. Place the earpiece securely in your left ear.



Tip

Start with a low intensity and increase gradually to a comfortable level.

Quick Start

- 1** Charge the device fully before first use.



- 2** Connect the earpiece to the device.



- 3** Lightly moisten the contact area.




- 4** Place the earpiece in the left ear and power on the device.



Starting a session

- Choose a mode that suits your needs.
- Press Start/Pause to begin the session.
- Increase intensity slowly to a comfortable level.

Earpiece Placement

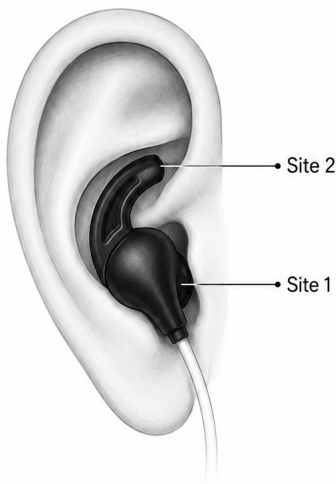

Placement guide

- Use the left ear.
- Ensure the earpiece sits comfortably and makes even contact.
- Adjust gently until the fit feels stable.
- Stop use if you feel significant discomfort.

Fit check

- ✓ Secure fit
- ✓ Even contact
- ✓ Comfortable intensity

Neuvago Wellness Device Pack



1. Select a mode

- Press the M button to cycle through the available modes (Sleep, Meditation, Relax, or Relief).
- The current mode icon will be shown on the display.

2. Start and adjust intensity

- Press the ▶/|| button to start or resume the session.
- Use the + and - buttons to adjust the intensity to your comfort level.
- Adjust gradually until the sensation is noticeable but comfortable.

3. Pause or end the session

- Press the ▶/|| button to pause. Press again to resume.
- Press and hold the power button for 2 seconds to end the session and turn off the device.



Comfort

A light tingling sensation may be expected.
Increase intensity only to a comfortable level.

Suggested routine

- ✓ Use in a calm seated position.
- ✓ Follow your guided routine.
- ✓ Return to regular sessions over time.

Care, Cleaning, and Storage

Cleaning

- Wipe the device surface with a soft, dry or slightly damp cloth. Do not use abrasive cleaners or solvents.
- Gently clean the earpiece with a soft, dry cotton swab or lint-free cloth after each use.
- Ensure the device is completely dry before storing or using again.



Storage

- Store the device in a dry, cool place away from direct sunlight
- Keep the device in the included storage pouch when not in use.
- Avoid placing heavy objects on top of the device or accessories

Important safety notes

- Do not immerse the device in water or any liquid.
- Do not use the device on broken, irritated, or damaged skin.
- Keep the device out of reach of children.
- Stop use and consult a healthcare professional if irritation or discomfort occurs.
- Use the device only as directed in this user manual.

Specifications & Support

Model:	NVG-001
Input:	5V ---1A
Battery capacity:	2000 mAh
Intensity:	0-8mA
Pulse width:	200-300 μ s
Frequency:	1-100 Hz

Support

Neuvago

www.neuvago.com

support@neuvago.com

Made in China

Compliance



Do not dispose of this device with household waste.